

Entrix NX



Pump

The Intelligent Pressure Sensing Technology (IPS), responds to patient movements on the mattress by automatically adjusting internal mattress pressure, this allows Entrix NX to regulate airflow and continuously provide total envelopment for our patients.



Alternating
Pressure Relief



Multiple Weight Settings



Whisper Quiet




CLP Therapy
Pressure Redistribution



Lightweight Power Unit



Memory Recall

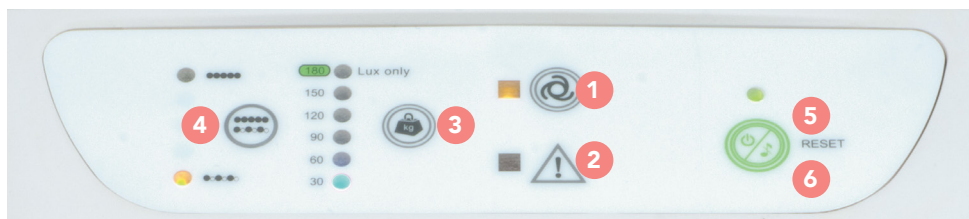

GENERAL FEATURES

1 Intelligent Pressure Sensing	Responds to patient movements on mattress by automatically adjusting internal mattress pressure.
2 Low Pressure Alarm	Audible and visible alert
3 Mattress Size Selection	Compatible with multiple mattress sizes which can be selected at start up
4 Dynamic Mode	2:1 alternating therapy, 10 minute cycle
4 Static Therapy/CLP	All air cells maintain constant pressure support
5 Power Failure Alarm	Audible and visual power failure alarms to alert caregivers in case of abnormal power outage
5 Memory Recall	Restores to the previous systems automatically after abnormal power outage
6 Alarm Mute	Disable the audible indicator system.
Transport Mode	Connect the two quick couplings together to retain air in the mattress, can stay inflated for up to 15 hours

TECHNICAL DETAILS

Risk/Category	High to Very High
Waterlow Score	15+ to 20+
Pressure Sore Grading	2 to 4
Safe Working Load	30 to 180kg
Compatible Support Surfaces	5" Overlay, 7" and Lux Mattress Replacements - Single and King Single
Travel/Weight	Pump weighs 2.3kg
Clinical Settings	Suitable for Hospital, Home Care and Nursing Home environments
Warranty	2 year

All Carilex Medical ANZ products are registered with the Therapeutic Goods Administration.